



Ampule Case (3 EA)



Main Body



Charging Cradle



Charging Adapter



User Manual

Components

Parts	Category	Unit	Specification	Remark
Main Body (+ Ample Case)	Weight	g	85	
	Injection Volume	ml	5 Stages (Checking with LCD on Main Body) 1/16 : 0.11 ml, 1/8 : 0.22 ml, 1/4 : 0.45 ml, 1/2 : 0.90 ml, FULL : 1.8 ml	
	Injection Time	s	Manual - 1st stage : 320s, 2nd stage : 160s 3rd stage : 30s Auto - 1st stage : 320s, 2nd stage : 170s 3rd stage : 48s	When injecting 1.8 ml
	Injection Mode	—	AUTO, MANUAL, USER	
	Battery	V	DC 3.7	Lithium Polymer
		mAh	550	
	Charging Time	h	6	
Charging Cradle	Size	mm	202.5 X 24.7 X 28	
	Weight	g	140	
Charging Adapter	Rated Input	V	AC 100 - 240	
		Hz	50 - 60	
	Rated output	V	DC 5	
		A	2.0	
	Weight	g	104	
Dental Needle(Disposal) *Separate Purchase	Size	mm	27G/16, 25, 30 30G/ 12, 16, 25 31G/ 12	

DENOPS

DENTIS No Pain System

DENOPS, the pain-free
anesthesia your patients
will enjoy!



CHANGE your
CLINIC



Tel. +82-1899-2804 | Fax. +82-53-583-2806
99, Seongseoseo-Ro, Dalseo-Gu, Daegu, Korea



Tel. +1-323-677-4363 | Fax. +1-323-677-4366
6 Centerpointe Drive, Suite 600 La Palma CA 90623



Tel. +021-5111-3828 | Fax. +021-5111-3828
上海市长宁区中山西路933号2205室

©2021 DENTIS Co., Ltd.
All Right Reserved. Specifications are subject to change without notice.
Trademarks are the property of DENTIS Co., Ltd. Or their respective owners.
D-L-DE-V1-202109-ENG

DENOPS

DENTIS No Pain System

- Light and slim design for easy use
- Injection speed control (Low/Med/High)
- User Mode makes it even easier to operate
- Voice assistance


Remove air bubbles

Injection and aspiration

LCD menu selection and mode change

LCD

- Battery level, injection volume and speed display
- Checking injection mode

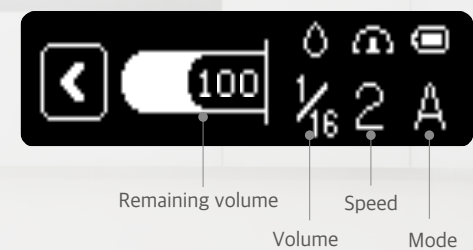

Wireless charging using a dedicated charging cradle

- Inclined design prevents leakage of anesthetic solution
- Battery level indicator



Plastic Ampule Case With Silicon Control Button

- Controlling injection speed and suction function
- Visual checking of injection volume and suction function with transparent Ampule case and window
- Autoclave available material for preventing cross infection